Tendri

PACK DUO MÉDICAL

Frontal sans contact & numérique

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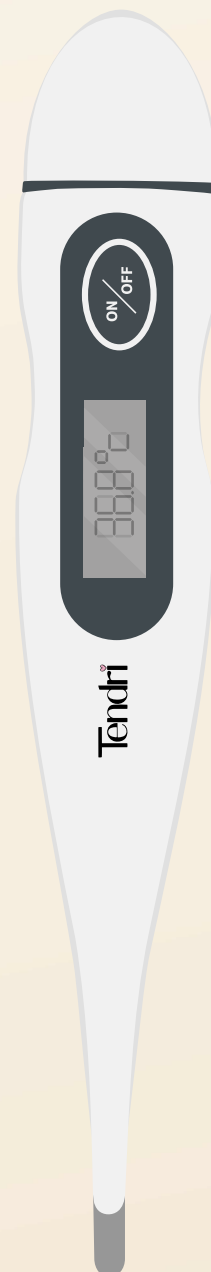
1

Power on

Press the Measure button to switch on.

2

Select age

👶 Newborn (0–3 months) 👶 Baby (3–36 months) ⓘ
Adult (+36 months)

3

Position

Aim at the centre of the forehead, 0–3 cm away, perpendicular to the skin.

4

Measure

Briefly press the Measure button.

5

Read result

The temperature appears instantly, in under a second.

Fever Alert System

GREEN

Normal temperature — No concern

ORANGE

Moderate fever — Monitoring recommended

RED

High fever — Consult a health professional


Memory — Short press on MEM: recall the last 35 readings.

°C/°F unit — Hold AGE for 3 sec to switch.

Silent mode — While on, hold MEM 3 sec to mute/unmute.

Object mode — When off, hold MEM 6 sec: bottle, bath, room (0–100 °C).

Fine tuning — Hold MEM + AGE 5 sec: ±2 °C, on medical advice.

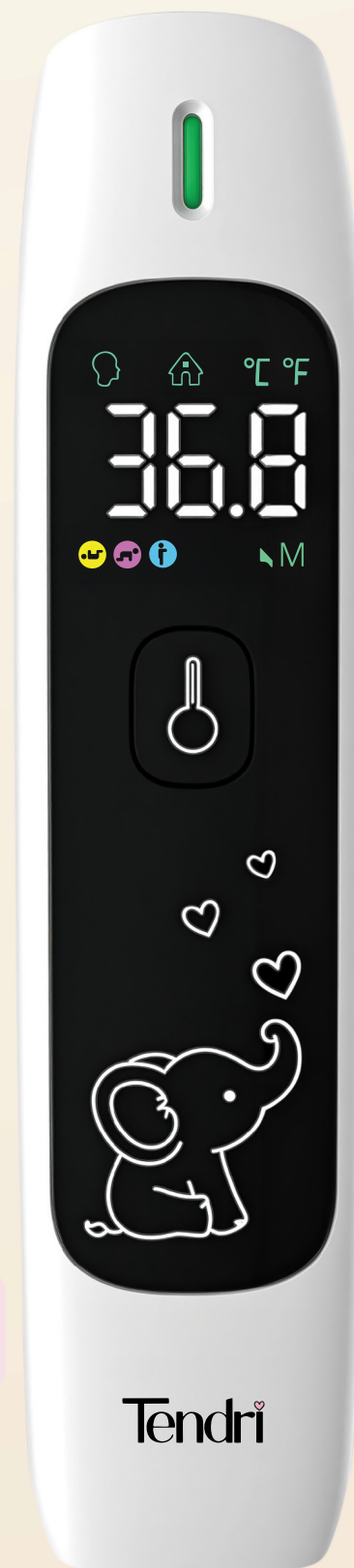Clean the forehead with a clean tissue before measuring to remove sweat or impurities.

Wait 30 min after physical activity, a bath or a meal for a reliable reading.

♥ Digital Thermometer

1

Power on

Press the ON/OFF button. A beep sounds, the screen lights up.

2

Choose the method

Oral: under the tongue • Underarm: under the arm • Rectal: use the probe (baby < 3 months).

3

Measure

Keep the thermometer in place until the beep.

4

Read result

The temperature shows on screen.

5

Turn off & clean

The thermometer switches off automatically after 8 min.
Clean the probe with a soft cloth and water.

i

Tip: clean the probe before and after each use. Do not immerse the thermometer in water.

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EN

- Use only as directed.
- Do not disassemble or modify the device.
- Keep out of reach of children.
- Batteries can be dangerous if swallowed.
- If in doubt, consult a doctor.
- Do not expose to heat, humidity or direct sunlight.

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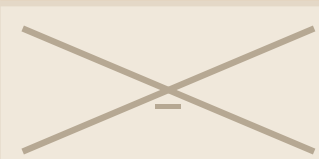
Infant mode: no orange level — goes straight to red from 37.5 °C.

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For guidance only — not a substitute for medical advice.

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Baby under 3 months with high temperature: seek advice promptly.

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EN

If the problem persists, contact customer service.

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ENGLISH

*Caring with
tenderness, measuring
with precision*

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
All rights reserved © Tendri

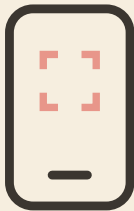Club Tendri

Inclus • Included • Inklusive • Incluso • Incluido

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EN Scan to record temperatures and medication.
Supplementary tool, no diagnosis.

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